

OSOM[®]

Influenza A&B Test



FOR USE WITH NASAL SWAB SPECIMENS


FEATURES:



- ✓ **ACCURATE, DIFFERENTIAL DIAGNOSIS;** Differentiates between influenza A&B
- ✓ **QC INSIDE;** Includes 2 additional test sticks for external QC testing
- ✓ **NASAL SWAB SPECIMENS**
- ✓ **FAST;** Less than 1 minute hands on time; results in 10 minutes or less
- ✓ **ROOM TEMPERATURE STORAGE**




TEST PROCEDURE / RESULTS

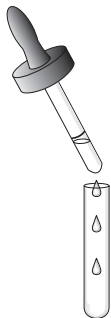

Add extraction buffer to the tube.



10x
Add swab to tube and vigorously mix.

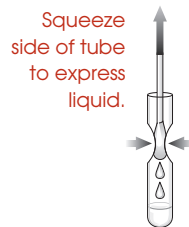


Squeeze side of tube to express liquid.

Place the test stick into test tube.



Read results at 10 MINUTES.

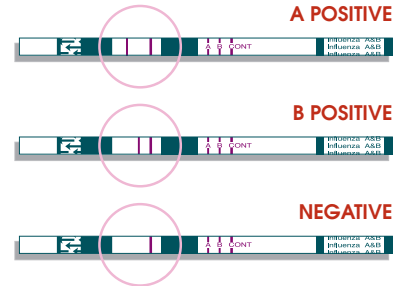

ORDERING INFORMATION

○ OSOM[®] INFLUENZA A&B TEST
INFLUENZA A&B CONTROL KIT

CONFIGURATION

25 TESTS + 2 TESTS FOR QC

CATALOG NUMBER

190 OR 190E (CE MARKED)
191 OR 191E (CE MARKED)



THE AMERICAS
Sekisui Diagnostics, LLC
31 New York Avenue
Framingham, MA 01701
Phone: 800 332 1042
Fax: 800 762 6311
Email: questions@sekisuidiagnostics.com
ISO 13485:2003 Certified

INTERNATIONAL
Sekisui Diagnostics (UK) Limited
50 Gibson Drive
Kings Hill, West Malling, Kent, ME19 4AF, UK
Phone: +44 (0) 1732 220022
Fax: +44 (0) 1732 220024
Email: info@sekisuidiagnostics.com
ISO 9001:2000 Certified, ISO 13485:2003 Certified



OSOM[®], the OSOM Training logo and the QC Inside logo are registered U.S. trademarks of Sekisui Diagnostics, LLC. ©2011 Sekisui Diagnostics, LLC – All rights reserved.